



Introducing the

KaVo MASTERtorque™ Mini

The Best Gets Better



**11%
MORE
TORQUE**

Superlative Performance

MORE POWER

19 watts peak power (2,8 bar) for approx. 11% more torque* and stable cutting over different drive pressure.

PERFECT ACCESSIBILITY

The Mini Turbine runs with miniature burs and standard burs.

REDUCED RUNNING NOISE

Effective frequency optimization ensures a pleasant turbine noise and low noise levels of only 59 dB(A) with better acoustic.

*Compared to its KaVo predecessors

**On average

**NOW WITH
DIRECT
STOP
TECHNOLOGY**

Cutting Edge Technology

SAFE PREPARATIONS

Innovative Direct Stop Technology stops the bur in less than one second,** ensuring increased safety during treatments.

WORK FASTER AND MORE EFFICIENTLY

DST stops the bur in less than one second and contributes to the reduced drawback of oral cavity fluids and dusts.

SAFE HYGIENE


Shortest stopping time prevents retractive suction of aerosols and particles into the ball bearings for more hygiene and overall handpiece longevity.

Bur stopping time



Stopping time compared to other turbines

Retractive suction from the head

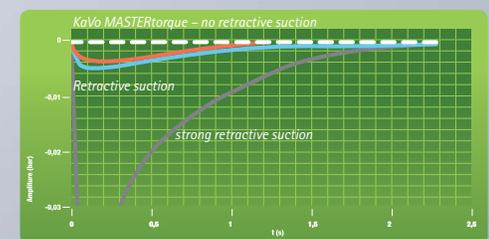

Go to **TRYKAVO.COM** to start your free and easy trial today!




KaVo. Dental Excellence.

888-ASK-KAVO · www.kavousa.com

Superior Ergonomics:

BETTER VISIBILITY

- Reduced head diameter on the bur side. Head height of only 12.1 mm
- Even with a 19 mm bur, the total head height is only 20.7 mm
- With KaVo's patented 100°/19° head and knee angle, molar access and visibility has never been easier

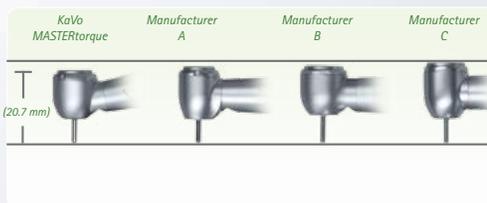

MORE LIGHT

- Preparation areas optimally illuminated at all times
- Preparation areas optimally visible in spray mode, light and spray nozzles are offset



MASTERtorque™ and MULTIflex™ Coupler Specifications

Product	Description	Part No.
MASTERtorque M8900 L	Standard head, 23W, Direct Stop Technology, ceramic bearings, 57 dB, 4-port spray, MULTIflex LUX connection	1008.7800
MASTERtorque Mini M8700 L	Mini head, 19W, Direct Stop Technology, ceramic bearings, 59 dB, 3-port spray, MULTIflex LUX connection	3001.0000
KaVo MULTI LED Bulb	MULTI LED bulb for 6-pin couplers	1007.5372
MULTIflex 460LED	6-pin coupler with KaVo MULTI LED	1007.3201
MULTIflex 465LED	6-pin coupler with KaVo MULTI LED bulb and spray regulator	1007.3197
MULTIflex 1390	5-hole coupler	0553.1390
QUATTROcare™ Plus	Automatic handpiece maintenance and purging system	1008.3805



KaVo QUATTROcare Plus

Prior to every sterilization, the unique cleaning and maintenance mixture produces highly effective cleaning foam inside the instrument. The automatic dispenser delivers the precise amount of cleaner and lubricant time after time, leading to longer handpiece life and less repair needs.

KaVo Dental
11727 Fruehauf Drive
Charlotte, NC 28273
888.ASK.KAVO (888.275.5286)
Fax 1.704.927.0816
www.kavousa.com

©2015 KaVo Dental Corporation. 1.011.3898/rev00/01.15



KaVo. Dental Excellence.

888-ASK-KAVO · www.kavousa.com